

ELECTRICAL OPERATING TABLE

ET300C

ET300C can meet various requirements from all department surgery. Extra wide tabletop, Long horizontal sliding that can be suitable for both X-ray and C-arm use. Adopted micro touch remote control that enables flexible and smooth adjustments on head plate, back plate and seat plate. with built-in kidney bridge. High Automation, low noise, high reliability. key parts adopted imported one, can be regarded as an Ideal electric operating table



Model \ Site	Length	Width	Tabletop lowest height	Elevation range
ET300C	≥2070mm	550mm	≤700mm	300mm

PARAMETER ET300C Electric Operating Table

Movements		Model:ET300C
Tabletop trendelenburg	Reverse Trendelenburg	≥ 25°
	Trendelenburg	≥ 25°
Lateral Tilt (Left And Right)		≥ 15°
Head plate	Up	≥ 45°
	Down	≥ 90°
Back plate	Up	≥ 75°
	Down	≥ 20°
Leg plate	Up	≥ 15°
	Down	≥ 90°
	Outward	≥ 90°
Kidney Bridge Up		≥ 120mm
Horizontal sliding		≥ 300mm

Product Description

Engineered for increasingly complex surgical scenes, ET300C comes with intelligent automatic control systemsatisfying various demands of general surgeries. Extra wide tabletop, extreme loading system,long horizontalsliding that can be suitable for both X-ray and C-arm. ET300C is anideal electric operating table.

- I. Longitudinal slide up to 300mm for free access to C arm and X-ray Transparent
- II. Extreme weight load capacity of 275kgs with highest safety and stability
- III. Memory pad with antistatic ,waterproof design
- IV. Standard override panel controlmakes the control safer
- V. Battery power for 50-80 operations



- ✓ Bed Frame 1pc
- ✓ Anesthesia Screen Frame 1pc
- ✓ Mattress 1set
- ✓ Kidney Bridge Handle 1pc
- ✓ Arm Support with Cushion 2pcs
- ✓ Leg Support with Cushion 2pcs
- ✓ Shoulder Support with Cushion 2pcs
- ✓ Body Support with Cushion 2pcs
- ✓ Foot Support with Cushion 2pcs

Leg Support with Cushion



Arm Support with Cushion



Body Support with Cushion



Shoulder Support with Cushion



 READY	 UP	 DOWN	 Bluetooth connect fail
Standby	Tabletop up	Tabletop down	Bluetooth connect fail


 TREND	 REV.TREND	 TILT L.	 TILT R.
Trendelenburg	Reverse	Lateral tilt left	Lateral tilt right

 BACK UP	 BACK DOWN	 SLIDE	 SLIDE
Back plate up	Back plate down	Horizontal sliding	Horizontal sliding



LEG UP

Leg plate up together

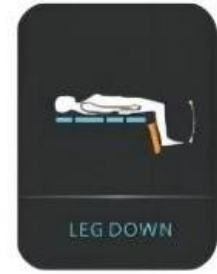

LEG DOWN

Leg plate down



FLOOR LOCK

Brake lock



FLOOR UNLOCK

Brake unlock


Powder Supply

Internal battery maintaining power supply for 50-80 operations and with low battery warning function.



LEG UP

Left electric leg up

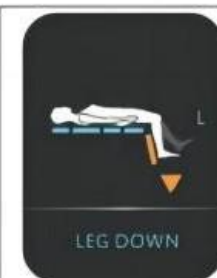


LEG DOWN

Left electric leg down



LEG UP

Right electric leg up

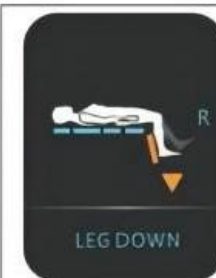

LEG DOWN

Right electric leg



Head Plate

Detachable double head plate and leg plate meet operation position needs and use Germany cylinder assist system.



LEVEL

Zero position



REFLEX

Reflex position



FLEX

Flex position




CHAIR

Chair position

Wash Basin Hole

Five section table top including Head, Back, Shoulder, Hip and Leg, all section applicable for X ray use and under each section with X ray cassette layer.

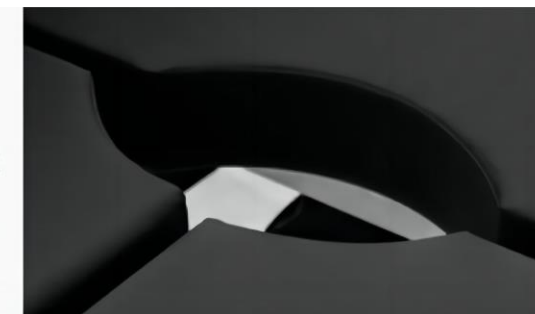

Menu



Custom setup



Maintenance



Battery information



Kidney Bridge

Built-In Kidney Bridge, More Convenient And Nice Appearance.



Remote Panel

All function controlled by hand controller in the column, and can be locked.

C-arm

Extra-Long Horizontal Sliding, Available For C-Arm Operation.

